



PERIMETER
PERISMART

PLUG & TOUCH BY RODENSTOCK

Compact design.
Portable and lightweight.



- High reproducibility
- LED technology
- Easy and fast operation
- Early Glaucoma detection
- Multiple illustration possibilities
- Small, compact and portable
- Reliable fixation control / digital eye tracking


QUALITY IN DETAIL

A literally 'easy-to-handle' solution for modern perimetry and screening – the Perismart is a compact, portable and mobile device combined with the advantage of a hemispherical perimeter. Shielded insight, electronic controlled chin rest and last but not least a structured and user optimised application software are only a few out of many other nameable points.

The software operation plus features as customised menu configuration and follow-up pattern appliance make the Perismart an ideal working equipment.

Test fields

The forthcoming test fields allow an extensive Glaucoma screening besides the standard testing.

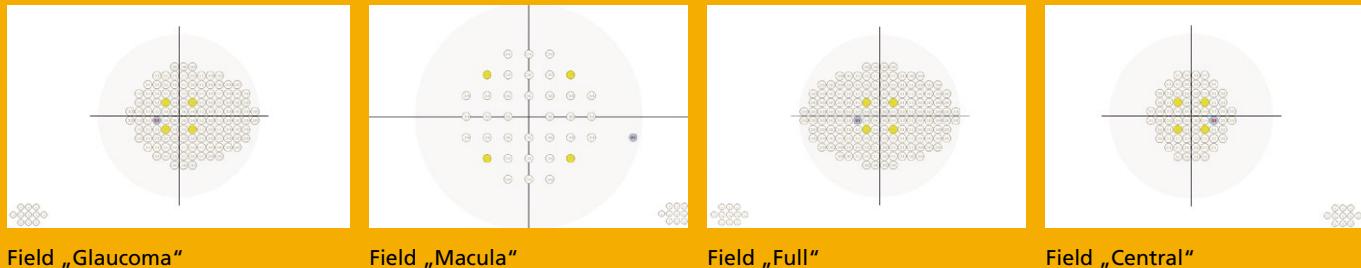

Field „Glaucoma“


Field „Macula“

Field „Full“


Field „Central“

Occupational medicine

The Perismart complies with the demands of industrial and occupational medicine (FeV) and shows test fields which correlate according to these requirements.

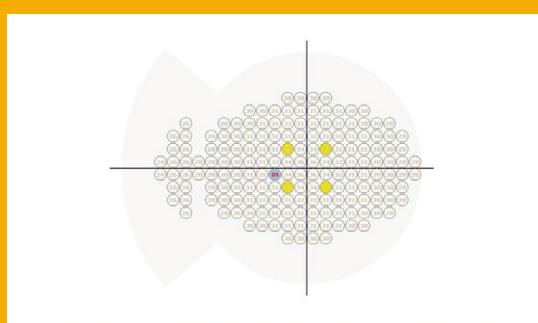

Field occupational medicine

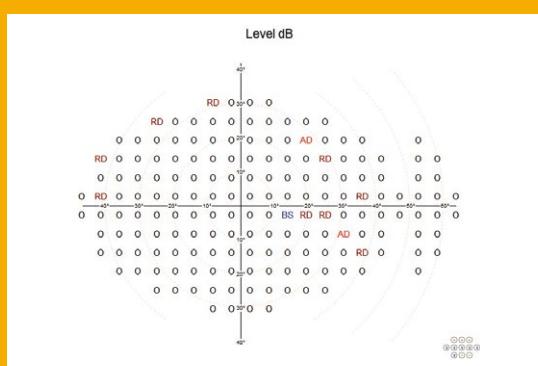

Result screen occupational medicine

SPECIFICATIONS

Stimulus bowl type	Part hemispherical
Bowl radius	300 mm
Eccentricity	50° horizontal (+30° mov. fix. point) 30° vertical (+10° mov. fix. point)
Background illumination	10 cd/m ² (32asb) auto level control
Stimulus	LED, white, Goldmann III
Stimulus intensity	0.01 cd/m ² (0.03 asb) to 333,33 cd/m ² (1000 asb) 3dB/1dB
Exposure time	0.1 to 9.9 sec adjustable
Chin rest adjustment	Electrically controlled
Fixation control	Heijl Krakau (blind spot method) and CCD camera (digital eye tracking)
User interface	Multilingual

Dimensions & weight

WDH	548 x 450 x 382 mm
Weight	Approx. 9 kg

Electrical requirements

Input	115 to 230 V
Frequency	50/60 Hz
Power consumption	Max. 60 W

